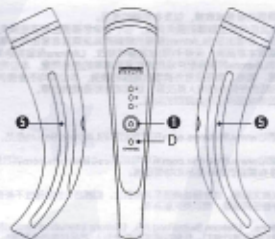
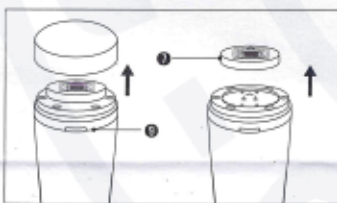
RF-300

A

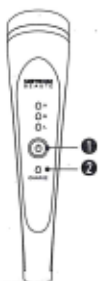
B

C

D

Thank you for purchasing this Lifetrans® product. This user guide outlines the operating procedures and safety precautions. Please read it carefully to make sure you are using your product properly and retain it for future reference.

How to Use

A. Radio Frequency & Red Photon Light Treatment

1. Remove any makeup before use. Apply skincare cream (not included) onto the treatment area. For best results and conductivity, you may choose to use a special radio frequency gel or cream (not included). Press the power button ① 1 time to turn on the device. The default intensity level will start with L (low).
2. Hold the device with both sides of your hand touching the conductive pole side strips ④ to activate the radio frequency mode. The radio frequency massage head ⑤ must always be in full contact with the skin to activate the red photon light ⑥ and start the treatment. Slowly and gently glide the device over your entire face with light pressure. You can use the device to target fine lines and wrinkles while the red photon light works simultaneously to promote collagen and elastin. You may choose your radio frequency treatment method: the default level L (low), M (medium) or H (high) depending on your skin's needs and sensitivity.
3. The device will automatically shut down after 15 minutes of continuous use. If you want to shorten the usage time, simply press and hold ① for 3 seconds and the device will shut down. After use, gently clean the massage head with a cotton pad.

Power & Intensity Level Settings

1. Press the power button ① 1 time: Device turns on and defaults to intensity level L (low)
2. Press ① 2 times: Switch to intensity level M (medium)
3. Press ① 3 times: Switch to intensity level H (high)

Note: The red photon light will activate whenever the device is turned on and makes contact with your skin.

Note: You may have a warm, stinging or pricking sensation during the treatment, depending on different levels of skin sensitivity. Some users also have mild skin redness after the treatment. These effects usually resolve after a short amount of time. If your skin is unusually hot during the treatment, please stop treating immediately and let your skin cool. Always keep your skin hydrated with moisturizers after use.

B. Cleaning the Device

1. Do not use any kind of liquid to clean the product. Keep in mind that the surface of the product should only be wiped with a clean dry cloth.
2. It is recommended to release the radio frequency massage head to clean after use. Press the release button ② at the back of the device to release the massage head ⑤.

C. Recharging the Device

1. Plug the included DC charging cable ① into the charging port ④ on the device.
2. Connect the other end of the DC charging cable ② to a power supply or wall socket to begin the charging process.

Note: Do not move back and forth in a push and pull motion. Move the device in one direction only. For best results, use up to 4 times a week for 3 minutes per session.

D. LED Indicator

Low battery: Blue LED light flashes once per second (only visible when power is on) ①

When charging: Solid red LED light appears ②

Battery fully charged: Solid red LED light disappears ③

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- Safety Precautions**
Please follow the rules outlined in this guide closely to ensure your safety and flawless operation:
- This product is only intended for use on the face.
 - Discontinue using the product immediately if any abnormal situation, allergic reaction or discomfort occurs during use.
 - This product is a skin care device and is not made for medical purposes. If you are under medical treatment, please consult with your doctor before use.
 - This product is not suitable to use on body parts that are embedded with silicone, metal or other foreign objects. Do not use directly on the eyes.
 - It is recommended to wait at least 6 months after any facial treatments such as a face lift or Botox injections before using the product. If you are unsure, consult with your doctor.
 - Thoroughly clean the product after use.
 - Do not use any kind of liquid to clean the product. Keep in mind that the surface of the product should only be wiped with a clean dry cloth.
 - This product has a special lacquer - when cleaning, do not use any abrasive solution which may cause damage to the lacquer or appearance.
 - Repair, removal or modification of this product is strictly prohibited. Improper handling may result in short circuits, product damage or even fire. Maintenance and warranty services will be void if unauthorized changes have been made to the product. Please contact an authorized dealer to obtain after-sales service if the product has a problem.
 - Be cautious of any edges, uneven surfaces, metal parts, accessories or packaging to prevent possible injury or damage. Take note that this product should not be placed in your mouth.
 - Do not use this product in locations with high temperatures or high humidity levels (e.g., bathrooms). Ideally, the product should be used in a location free of excessive dust and not exposed to direct sunlight.
 - Do not insert foreign objects or liquids into the product, this could damage the product.
 - Refrain from dropping the product or bumping it against other objects. Avoid damage, collision and excessive pressure by always handling the product with care.
 - Do not place this product near strong magnetic or electric fields - this may cause product failure to occur.
 - Do not use this product in hazardous locations.
 - The product may explode if disposed of in a fire.
 - Remove the product from the power supply within 5 hours of the finished charging cycle.
 - Do not use internally and avoid using near center of the throat to avoid thyroid stimulation.
 - Do not use on open sores, wounds or infected areas.
 - Never use internally (e.g. in mouth, nose, etc.)
 - Do not use if pregnant, have a pacemaker, heart disease history or other medical condition.

Children
Keep out of reach of children. Do not allow children to play with your product or its accessories. They could hurt themselves or others and could accidentally damage the product. Your product and its accessories may contain small parts that could become detached and create a choking hazard.

Disposal of Electrical and Electronic Equipment
This symbol indicates that all electrical and electronic equipment included with this product should not be treated as household waste. Instead, it should be left at the appropriate collection point for the recycling of electronic equipment, by ensuring this product is disposed of correctly, you are helping to prevent negative consequences for the environment and human health. To recycle this product, please contact your local city office, environmental protection agency, household waste disposal service or where you purchased the product to learn more about appropriate recycling methods.

Battery
It is recommended that you completely use the initial charge up, then recharge the product at least two more times to ensure the best possible performance.

A rechargeable battery can be charged and discharged over one hundred times. However, it will eventually wear out - this is not a defect and corresponds to normal everyday use. When the operating or standby time is noticeably shorter, it may be a sign that the battery has reached the end of its life span. Please do not make repairs, modifications, or installations without contacting an authorized Lifetrons representative.

Limited Warranty
If within 12 months of purchase this product or any part thereof is proven to be defective by reason of faulty workmanship or defect in materials, we will repair or replace the same product free of charge of labor or materials on condition that:

- The original invoice or sales receipt with the date of purchase is presented with the defective unit.
- Lifetrons does not take any responsibility for lost information of any kind and will not reimburse you for any such loss.
- The limited warranty is only valid if the original proof of purchase for this product, issued by a Lifetrons authorized dealer, specifying the date of purchase, is presented with the product to be repaired or replaced. Lifetrons reserves the right to refuse warranty service if this information has been removed or changed after the original purchase of the product from the dealer.
- The warranty does not cover product failure which has been caused by the use of accessories or other peripheral devices which are not Lifetrons branded original accessories intended for use with the product. This warranty does not cover any damage caused by misuse, abuse, neglect or accidental damage.
- The warranty does not cover product failure caused by modifications, repairs or disassembly by any person not authorized by the manufacturer.
- Any defective unit or part shall become property of the manufacturer.

Additional Information
Register your Lifetrons Beauty product online at www.lifetrons.com.

Please visit our website www.lifetrons.com or contact our customer service center at cs@lifetrons.com to learn more about our new products and special offers.

This user guide was originally composed in English and was subsequently translated into other languages. The fidelity of subsequent translations cannot be guaranteed. In case of conflict between the English version and another language version, the English version takes precedence.

Product specifications are subject to change at the sole decision and discretion of Lifetrons Switzerland AG, Lifetrons International Ltd. or Li-Tek, without prior notice. Any devices or accessories not listed under "Official models" are in reference to compatibility only and are not included. Product images are for reference purposes only, color may vary from the actual product.

"Lifetrons Beauty", "Lifetrons Beauty Switzerland", "Lifetrons Beauty Professional Skincare", "Lifetrons Switzerland", "Beauty Switzerland", "WYS" are trademarks or registered trademarks owned by Lifetrons Switzerland AG, Lifetrons International Ltd. or Li-Tek. All other brands and product names listed herein are trademarks or registered trademarks owned by their respective holders and are hereby acknowledged and recognized as such. ©Copyright 2016 by Lifetrons Switzerland AG, Lifetrons International Ltd. and Li-Tek. All rights reserved.

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"Lifetrons BEAUTE", "Lifetrons BEAUTE SWITZERLAND", "Lifetrons BEAUTE PROFESSIONAL SKINCARE", "Lifetrons SWITZERLAND", "BEAUTE SWITZERLAND", "WYS" are trademarks or registered trademarks owned by Lifetrons Switzerland AG, Lifetrons International Limited and Li-Tek. All other brands and product names listed herein are trademarks or registered trademarks owned by their respective holders and are hereby acknowledged and recognized as such. ©Copyright 2016 by Lifetrons Switzerland AG, Lifetrons International Limited, and Li-Tek. All rights reserved.

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