


Ultimate Facial Kit

Skin Care Device with 6 Functional



RP-100

All pictures shown are for illustration purpose only, actual products may vary.
圖片僅供參考，實物為主。 圖片僅供參考，實物為主。

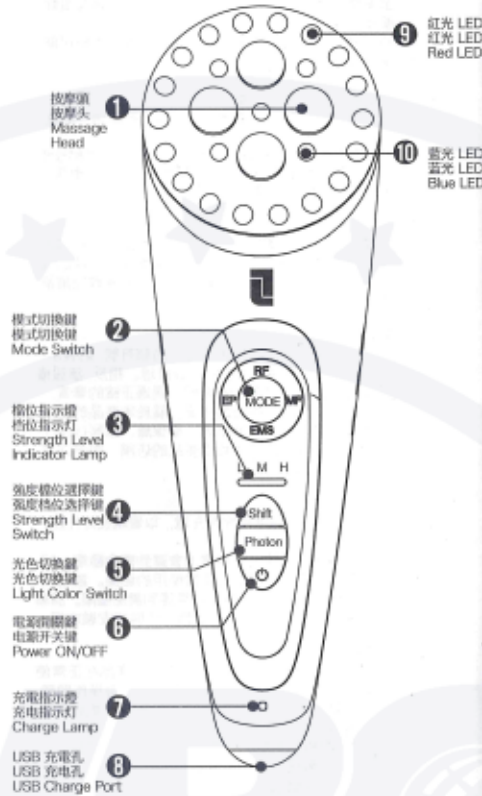

EN User Guide

Thank you for purchasing this Lifetrans® product. This user guide outlines the operating procedures and safety precautions. Please read it carefully to make sure you are using your product properly and retain it for future

How to Use

A. EMS + Vibration mode

1. Click ⑥ to start the device, light EMS and L light, EMS + vibration function is on standby mode, accompanied by a Di sound.
2. Press ④ to select the strength level, there are L for low, M for medium, H for high, three different levels can be selected, each with a Di sound.
3. Hold onto the device, place ① gently onto the skin and start treatment care.

B. EP mode

1. Press ⑥ to start the device, press ② and the indicator lights for EP and L will light up, EP function will be on, and every press will have a Di sound.
2. Press ④ to select the strength level, there are L for low, M for medium, H for high, three different levels can be selected, each with a Di sound.
3. Hold onto the device, place ① gently onto the skin and start treatment care.

C. RF Mode

1. Press ⑥ to activate the device, press ② twice and the indicator lights for RF and L will light up, RF function will be activated and every press will have a Di sound.
2. Press ④ to select the strength level, there are L for low, M for medium, H for high, three different levels can be selected, each with a Di sound.
3. Hold onto the device, place ① gently onto the skin and start treatment care.

D. MF mode

1. Press ⑥ to start the device. Press ② three times. When the indicator lights for MF and L light up, the MF function will be activated and each press will have a Di sound.
2. Press ④ to select the strength level, there are L for low, M for medium, H for high, three different levels can be selected, each with a Di sound.
3. Hold onto the device, place ① gently onto the skin and start treatment care.

E. Red. Blue Light mode

1. Press ⑥ to start the device. Press ③ to turn on the red LED ⑨. Press ⑤ and ⑩ for blue LED to light up, there will be a Di sound for every press.
2. Hold onto the device, place ① gently onto the skin and start treatment care.

F. Cleaning of device

After using the device, please use a dry cloth or paper towel to clean and store it.

G. How to shut down

When the pre-set timing for treatment is up, the device will shut down automatically. If you want to shorten the use of time, just press ③, the function will stop and Di three times.

H. How to charge

1. Turn on the USB plug on ⑧, insert one end of the charging cable into ⑧, and connect the other end to the external power supply (5V power supply, such as computer USB port and mobile phone charger).
2. When charging, ⑩ blue LED light flashes, when fully charge, ⑨ LED light will stop flashing and in stable mode.

EN Safety Precautions

Safety Precautions

Please follow the rules outlined in this guide closely to ensure your safety and flawless operation:

- This product is meant for facial skin.
- This product is for skin care only and is not intended for medical use. If you are receiving any medical treatment, please consult your doctor first before using.
- This device is not intended for use in body part where silicone, metal or other foreign materials are implanted, nor can it be applied directly to the eyeball.
- If you have a face lift or botox injection, please complete the course and use this device only after 6 months.
- After the treatment, keep the device clean.
- Do not use any detergent or liquid to clean the surface of this device. This product can only be cleaned with a dry cloth.
- This product has a unique paint protection, do not use any corrosive cleaning solution or it may cause damage to the surface.
- Do not attempt to repair, disassemble or modify this unit. Improper handling may result in short-circuiting of wires, product damage or even fire. Any unauthorized modification to the product, it will lose the function effectiveness and void of warranty. If the device is having problem, please go to the designated after-sales service centre.
- Any cutting edges, uneven surfaces, metal parts, accessories or packaging that may cause bleeding or injury by accident, be sure to use with caution. Do not expose the device to mouth.
- Do not use this device in places subject to high temperatures or high humidity, such as in bathrooms. Use this device in locations with low dust. Do not expose it to excessive heat such as sunshine, roset or similar.
- Do not put foreign objects or liquids in the product, which may damage the product.
- Do not put foreign objects or liquids in the device, which may damage the device.
- Do not place the device in a strong magnetic field or an electric field as it may cause the device malfunction.
- Do not use in a hazardous environment.
- This product may cause an explosion if left in an open flame.
- This device is extremely energy-intensive. Do not aim the nursing window at sensitive areas such as eyes or heart.
- To extend the product life time, it is recommended to do the charging- discharging every month for the power charger which is rarely used. If the power charger will not be used for a long time, it is better to reserve 50% battery power.

Children

Keep out of reach of children. Do not allow children to play with your product or its accessories. They could hurt themselves or others and could accidentally damage the product. Your product and its accessories may contain small parts that could become detached and create a choking hazard.

Disposal of Electrical and Electronic Equipment

⚡ This symbol indicates that all electrical and electronic equipment included with this product should not be treated as household waste. Instead, it should be left at the appropriate collection point for the recycling of electronic equipment. By ensuring this product is disposed of correctly, you are helping to prevent negative consequences for the environment and human health. To recycle this product, please contact your local city office, environmental protection agency, household waste disposal service or where you purchased the product to learn more about appropriate recycling methods.

Battery

It is recommended that you completely use the initial charge up, then recharge the product at least two more times to ensure the

best possible performance.

A rechargeable battery can be charged and discharged over one hundred times. However, it will eventually wear out—this is not a defect and corresponds to normal everyday use. When the operating or standby time is noticeably shorter, it may be a sign that the battery has reached the end of its life span. Please do not make repairs, modifications, or installations without contacting an authorized Lifetrons representative.

Limited Warranty

If within 12 months of purchase this product or any part thereof is proven to be defective by reason of faulty workmanship or defect in materials, we will repair or replace the same product free of charge of labor or materials on condition that:

1. The original invoice or sales receipt with the date of purchase is presented with the defective unit.
2. Lifetrons does not take any responsibility for lost information of any kind and will not reimburse you for any such loss.
3. The limited warranty is only valid if the original proof of purchase for this product, issued by a Lifetrons authorized dealer, specifying the date of purchase, is presented with the product to be repaired or replaced. Lifetrons reserves the right to refuse warranty service if this information has been removed or changed after the original purchase of the product from the dealer.
4. The warranty does not cover product failures which have been caused by use of accessories or other peripheral devices which are not Lifetrons branded original accessories intended for use with the product. The warranty does not cover any damage caused by misuse, abuse, neglect or accidental breakage.
5. The warranty does not cover product failures caused by modifications, repairs or disassembly by any person not authorized by the manufacturer.
6. Any defective unit or part shall become property of the manufacturer.

Additional Information

Register your Lifetrons Beauté product online at www.lifetrons.com.

Please visit our website www.lifetrons.com or contact our customer service center at cs@lifetrons.com to learn more about our new products and special offers.

This user guide was originally composed in English and was subsequently translated into other languages. The fidelity of subsequent translations cannot be guaranteed. In case of conflict between the English version and another language version, the English version takes precedence.

Product specifications are subject to change at the sole decision and discretion of Lifetrons Switzerland AG, Lifetrons International Ltd. or Li-Tek without prior notice. Any devices or accessories not listed under "What's Inside" are in reference to compatibility only and are not included. Product images are for reference purposes only, color may vary from the actual product.

"Lifetrons Beauté", "Lifetrons Beauté Switzerland", "Lifetrons Beauté Professional Skincare", "Lifetrons Switzerland", "Beauté Switzerland", "MyB" are trademarks or registered trademarks owned by Lifetrons Switzerland AG, Lifetrons International Ltd. or Li-Tek. All other brands and product names listed herein are trademarks or registered trademarks owned by their respective holders and are hereby acknowledged and recognized as such. ©Copyright 2018 by Lifetrons Switzerland AG, Lifetrons International Ltd. and Li-Tek. All rights reserved.

簡中 安全守則

安全守則

為保安全和操作正確，請先詳細閱讀及遵守以下的安全守