

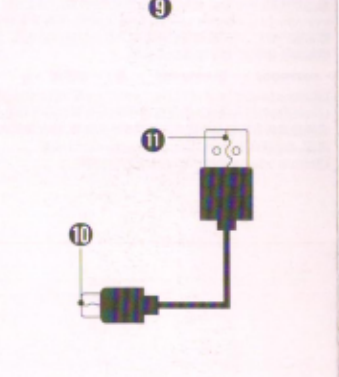
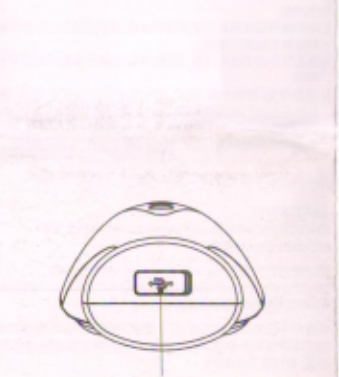
Beauty Lift V Face Device

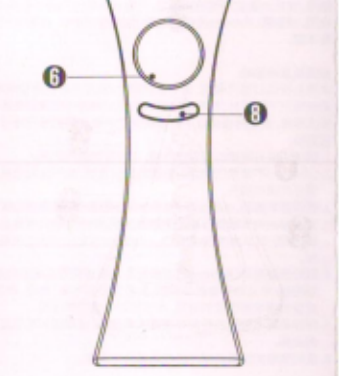TPML-100

(EN) All pictures shown are for illustration purpose only, actual products may vary.
(DE) Alle gezeigten Bilder dienen nur zur Illustration, die tatsächlichen Produkte können variieren.
(FR) Les images ne sont fournies qu'à titre d'information et le produit réel peut varier.
(ES) Las imágenes mostradas son sólo para referencia y el producto real puede variar.
(IT) Tutte le immagini mostrate sono a scopo puramente illustrativo, i prodotti reali possono variare.
(RU) Все изображения показаны только для иллюстрации, реальный продукт может отличаться.
(圖中) 圖片僅供參考，實物為主。

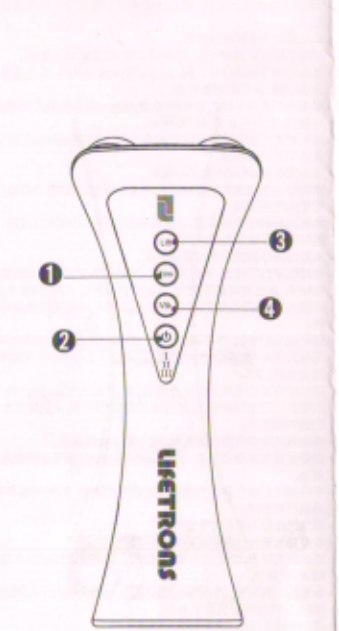Skincare Procedure

A: Clean

Clean your face with facial cleanser.

B: Lifting Treatment

Application Frequency: 15 minutes for each treatment, 2-3 times per week. It is recommended to use it with skin care products. When using the device, do not let the photon light shine directly on the eyes. The treatment technique is let the device close to the skin and treat according to the direction indicated by the arrow in the figure.

Massage from the chin to the ear back and then to the neck.

Massage your face upwards along the contours.

Massage from your eyebrows to the hairline.

C: Firming Treatment

Application Frequency: 15 minutes for each treatment, 2-3 times per week. It is recommended to use it with skin care products. When using the device, do not let the photon light shine directly on the eyes. The treatment technique is let the device close to the skin and treat according to the direction indicated by the arrow in the figure.

Massage outward and upward along the facial contour.

Massage from your eyebrows to the hairline.

Massage your eyes area in circle.

EN User Guide

Thank you for purchasing this Lifetrans® product. This user guide outlines the operating procedures and safety precautions. Please read it carefully to make sure you are using your product properly and retain it for future reference.

This product is applicable for face only. This device is only applicable to the face. Please do not let the photon light shine directly on the eyes when using it. It has the risk to cause damage to the eyesight.
Refer to instruction manual booklet

How to Use

Please remove makeup and clean your face before treatment. You could use the device with skin care cream or moisturizing cream.

- A. Power on / off**
 - Power on: press the power button ➊ for a short time, and its indicator will stay in blue light.
 - Power off: long press the power button ➋ for 2 seconds, and its indicator flashes twice in blue light.

B. Lifting treatment function(Warm,EMS & Red photon)

- Use the Lift button ➊
 - Press it for the first time to turn on the EMS function, and its indicator will stay in pink light. Press the second time to turn on the EMS and Warm functions, and the treatment head ➋ will start to warm. When the temperature reaches 40 °C, its indicator will flash in blue light, and you can start facial care.
 - Press it for the third time to turn on the Red Photon therapy function, and its indicator will stay in blue light. Press again to turn off the above functions. Max temperature 44±2 °C.
- Press ➋ for 3 seconds to open the EMS, Warm and Red Photon functions at the same time. After treatment, press ➋ to close it.
- Treatment time: Please follow the treatment process to operate the device. For facial treatment, suggesting using 15 minutes once. Please switch off the device if no usage.
- Micro vibration function can be used at the same time to increase the absorption of skin products.

- Press ➋ for 3 seconds to open the EMS, Warm and Red Photon functions at the same time. After treatment, press ➋ to close it.
- Treatment time: Please follow the treatment process to operate the device. For facial treatment, suggesting using 15 minutes once. Please switch off the device if no usage.
- Micro vibration function can be used at the same time to increase the absorption of skin products.

C. Firming treatment function(Comfortable temperature & Blue photon)

- Use the Firm button ➊
 - Press it for the first time to turn on the skin firming function, and the Firming massage head ➋ will start to work; when the temperature drops to 10 °C ± 2 °C, its indicator will flash in blue light, and you can start facial care.
 - Press it for the second time to turn on the Blue photon therapy function, and its indicator ➌ will stay in blue light. Press again to turn off the above functions.
- After the Firming function is turned on for 3 minutes, the Comfortable temperature function will be automatically turned off. Blue photon will be automatically turned off after turning on 15 minutes.

- Press Micro-vibration function ➊ once for continuous vibration mode.
- Press ➋ twice for intermittent vibration mode which simulated pat and massage.
- Press ➌ three times for strong and weak vibration mode which simulated push-pull massage.
- Press ➌ four times to turn off the vibration mode. In this cycle.
- It can be used alone or together with lifting or firming functions.

E. LED indicator

- Power button ➊ indicator: when working it stays lighting in blue, it will be closed when shutdown (under no charging state).
- Low power indicator ➋ will flash in pink light once every 2 seconds (only in the power on state).
- Full battery ➌ blue indicator light is on for a long time.

F. How to charge

- Insert the ➊ Micro-USB plug of the connected transformer (DC 5V,1A) into the Micro-USB port.
- When the device is charging, ➋ pink LED light will be on for a long time. After charging, ➌ blue LED light will be on for a long time. At this time, the device ➌ or ➋ can be disconnected.

G. Precautions

- Once the firming function is started, after completing the whole treatments processes, it needs to be shut down for 30 minutes to continue to use. If the interval is not long enough, it will affect the comfortable temperature.
- When charging, the device cannot be used. It can only be used when it is not charged.
- Red & blue Photon can't irradiate the eyes. When using the Photon function of the instrument, you need to wear the goggles.
- If no function is in use for 1 minute, the power indicator ➋ will flash twice in blue, and it will automatically shut down.
- The duration of each use of the firming function shall not exceed 3 minutes, and the function will be automatically turned off after 3 minutes. It's to protect the normal use of the electrical devices of the function from damage.
- In order not to affect the effect of the product, lifting and firming functions can not be used at the same time.

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